


## Force Scientific Technologies Pvt. Ltd.

### ANALYTICAL DEVELOPMENT LABORATORY

<b>Sequence Directory :-</b>	UHPLC3\Method Development\Sildenafil Tablet IP.seq		
<b>Sequence Name :-</b>	Sildenafil Tablet IP	<b>Instrument ID :-</b>	UHPLC3
<b>Injection Name :-</b>	Sildenafil Tablet IP	<b>Column No. :-</b>	n.a.
<b>Injection No. :-</b>	9	<b>Vial Position :-</b>	R:A7
<b>Injection Type :-</b>	Unknown	<b>Injection Volume :-</b>	20.0 µL
<b>Instrument Method :-</b>	Sildenafil Tablet IP	<b>Run Time :-</b>	20.0 mins
<b>Processing Method :-</b>	3D Quantitative	<b>Analyst :-</b>	adl
<b>Injection Acquired On :-</b>	08.04.2024 15:15:15	<b>Wavelength :-</b>	294 nm

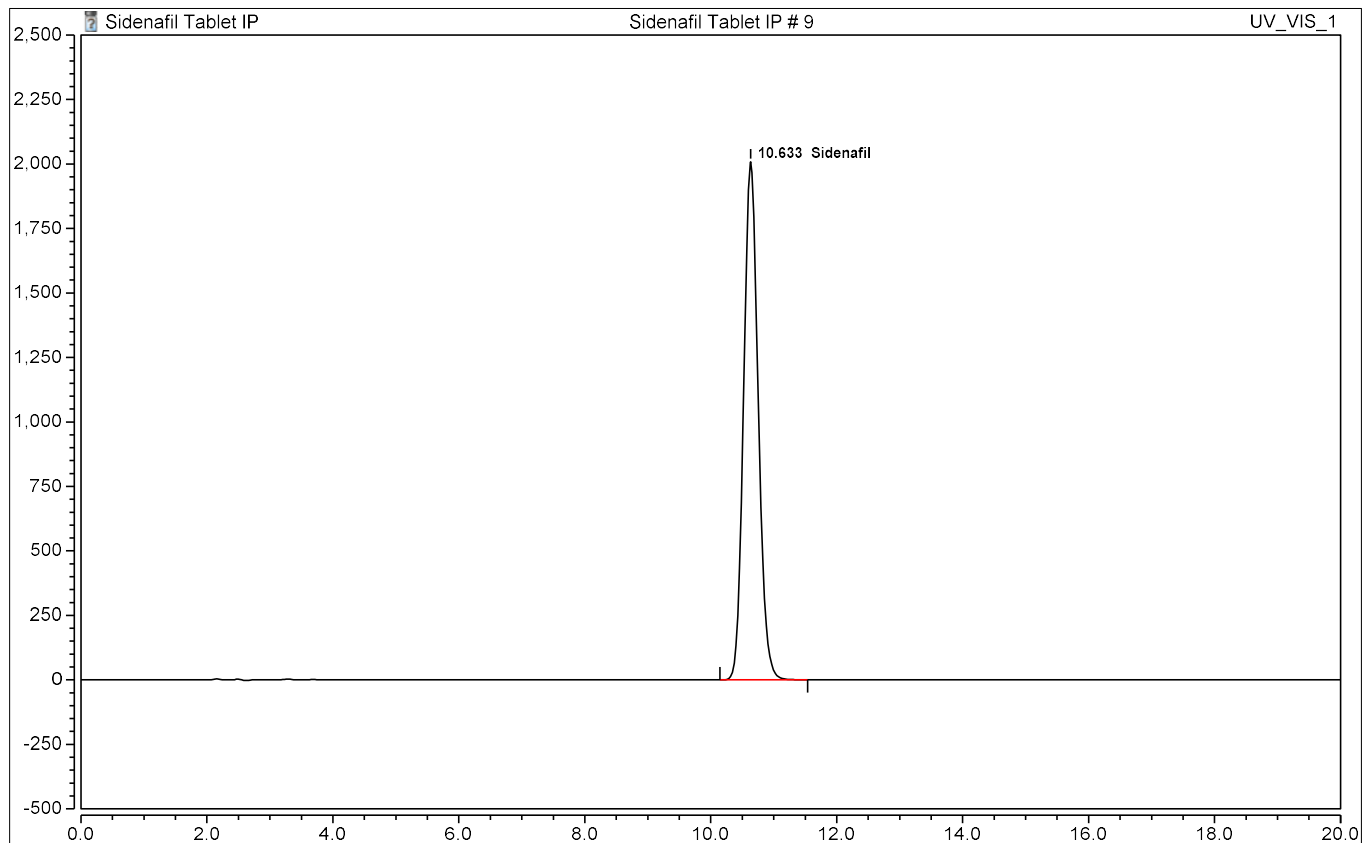

Peak No.	Ret.Time min	Peak Name	Area µAU*sec	Area %	Tailing (USP)	Resolution (USP)	Theoretical Plates	RRT
1	10.633	Sildenafil	32360693	100.00	1.11	n.a.	10192	n.a.
<b>Total</b>			<b>32360693</b>	<b>100.00</b>				